




Blood Glucose Testing Guide

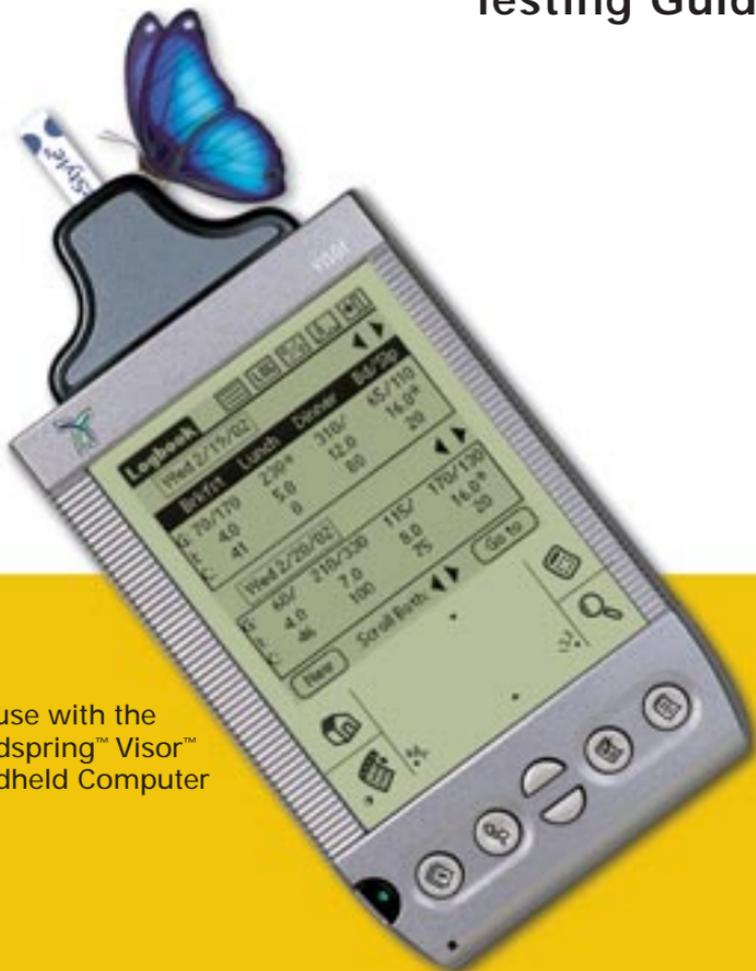

For use with the
Handspring™ Visor™
Handheld Computer

FreeStyle TRACKER™

Diabetes Management System

FreeStyle Tracker Diabetes Management System Blood Glucose Testing Guide

First, use the System Installation Guide to set up the FreeStyle Tracker System.

After the System is set up, use this Blood Glucose Testing Guide to test your blood glucose, or to perform a control solution test.

Detailed information about the material in this Guide can be found in the electronic FreeStyle Tracker System User's Guide on the CD-Rom.

How the FreeStyle Tracker System Should be Used	2
Important Health Related Information	4
Blood Glucose Test Overview	6
HandSpring™ Visor™ Handheld Computer (PDA)	7
FreeStyle Tracker™ System Contents	8
Step 1 – Prepare the FreeStyle™ Lancing Device	9
Step 2 – Insert the Module into the Visor PDA	12
Step 3 – Insert a Test Strip into the Module	13
Step 4 – Check the Test Strip Cal Code	14
Step 5 – Obtain a Blood Sample	16
Step 6 – Apply the Blood Sample	17
Step 7 – See the Test Result Displayed	19
Viewing Test Results in the Logbook	20
Viewing Test Results in the Diary	21
Logbook and Diary Screen Icons and Terms	22
Multiple Test Site Information	23
Additional Lancing Device Information	24
About FreeStyle Control Solution	25
Problem Solving Information	27
Error Messages	28
Caring for Your System	32
Specifications	32

The FreeStyle Tracker™ System is intended for use outside the body (i.e. *in vitro* diagnostic use only).

The FreeStyle Tracker System is for measuring blood glucose from samples taken from the forearm, upper arm, hand, thigh, calf, or fingers (for more information refer to the FreeStyle Tracker System User's Guide).

The FreeStyle Tracker System should only be used with FreeStyle™ Test Strips and FreeStyle™ Control Solution. Using other brands of test strips or control solution with the FreeStyle Tracker System may produce inaccurate results.

FreeStyle Tracker should only be used for:

- Testing blood glucose (sugar)
- Testing only with fresh capillary whole blood samples

FreeStyle Tracker should not be used for:

- The diagnosis of diabetes
- Testing on newborns
- Testing of arterial blood

WARNING

This System contains small parts that may be dangerous if swallowed. Keep away from young children.

If you are testing for hypoglycemia (low blood sugar), or if you suffer from hypoglycemia unawareness, we recommend that you test on your fingers.

CAUTION

Before using the FreeStyle Tracker System to test your blood glucose, please read all the instructions provided in this Blood Glucose Testing Guide and practice the testing procedures. Blood glucose monitoring should be done with the guidance of a healthcare professional.

Most healthcare professionals recommend testing when blood glucose values are stable, before meals and before bedtime. However, when blood glucose is changing, blood from the fingertips may show these changes sooner than blood from other sites.

IMPORTANT

The FreeStyle Tracker System is not designed to perform a blood glucose measurement while the system is placed in a Visor cradle or otherwise being charged. Remove the PDA from the HotSync Cradle or power source before doing a blood glucose measurement to avoid a system reset.

- Severe dehydration and excessive water loss may cause false low results. If you believe you are suffering from severe dehydration, consult your physician immediately.
- Test results below 60 mg/dL (3.3 mmol/L) mean low blood glucose (hypoglycemia). Test results greater than 240 mg/dL (13.3 mmol/L) mean high blood glucose (hyperglycemia).
- If you get results below 60 mg/dL or above 240 mg/dL, and you do not have symptoms of hypoglycemia or hyperglycemia, repeat the test. If you have symptoms or continue to get results that fall below 60 mg/dL or above 240 mg/dL, follow the treatment advice of your healthcare professional.
- If you are experiencing symptoms that are not consistent with your blood glucose test results AND you have followed all instructions described in this Blood Glucose Testing Guide, call your healthcare professional.

- Physiologic differences in the circulation between the finger and other test sites like the forearm, upper arm, thigh, calf and hand, may result in differences in blood glucose measurements from the other test sites and your finger. Differences in blood glucose measurements between the other test sites and your finger may be observed after eating, insulin medication or exercise.
- Changes in blood glucose may be observed in finger blood samples sooner than blood samples from the forearm and other alternate sites. Vigorous rubbing of the alternate test sites before lancing will help to minimize the difference between finger and alternate site test results.
- If you are testing for hypoglycemia (low blood sugar), or if you suffer from hypoglycemia unawareness, we recommend that you test on your fingers.

Times when you may be checking for low blood glucose include:


- When you are feeling symptoms such as weakness, sweating, nervousness, headache, or confusion.
- When you have delayed a meal after taking insulin.
- When your healthcare professional advises you to do so.

These are the main steps to test your blood glucose:

1. Prepare the FreeStyle™ Lancing Device, then set it aside.
2. Remove the cover from the upper rear of the Visor™ Handheld Computer (PDA) and insert the FreeStyle Tracker™ Glucose Measurement Module (if not already inserted).
3. Insert a FreeStyle™ Test Strip into the Module.
4. Check that the test strip Cal Code displayed on the Visor PDA matches the code on the test strip vial – if not, adjust the Cal Code on the Visor PDA display.
5. Clean and dry the site to be lanced, then obtain a blood sample using the FreeStyle Lancing Device.
6. Touch the test strip (that has been inserted into the Module) to the blood droplet until you hear a single beep or see the rotating “Clock” icon.
7. View the result on the display

The result is automatically saved in the Logbook and Diary. When the Logbook or Diary screens are showing, you can tap on the result and then tap the Details button to view more information, or to add, change, or delete information (such as comments). The FreeStyle Tracker™ System User’s Guide has more information about the Logbook and Diary.

NOTE: Words that are underlined (for example – Details) are on-screen commands or functions that you tap with the stylus to activate. Words that are capitalized have a special meaning or usage, or are product names.

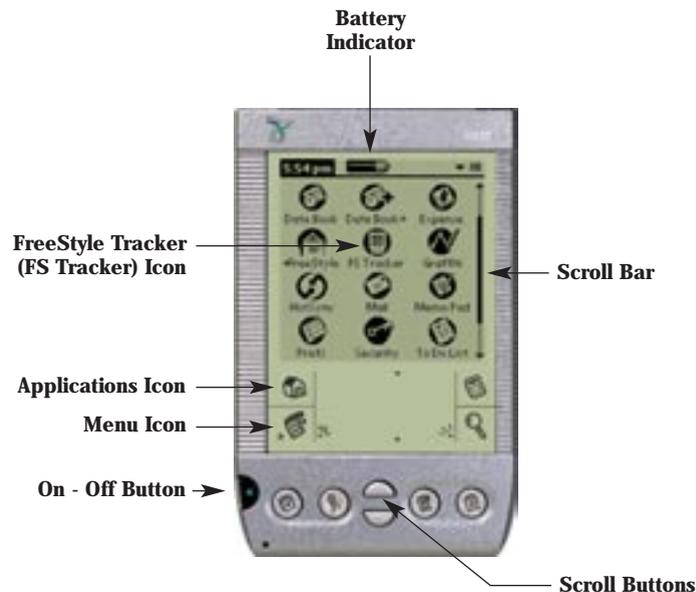

IMPORTANT

This Guide applies to FreeStyle Tracker Systems using the 1.2 program version only. To check which program version your System is using, while the FreeStyle Tracker System is running tap the  (Menu) icon, then tap About This Program. Do not use this Guide if it does not match the program version. Contact TheraSense Customer Service toll free at 1-888-522-5226 if you need a different version of this Guide.



* FreeStyle Tracker Software Kit Contents

** Not included in FreeStyle Tracker Monitoring Module and Software Kit

The FreeStyle Lancing Device



IMPORTANT

Lancing with the Clear Cap may be completed at the following test sites: forearm, upper arm, hand, thigh, calf.

Remove the Cap



1. Snap off the Clear Cap on the lancing device, like breaking a cracker in half.

Insert the Lancet



2. Insert a new FreeStyle™ Lancet firmly into the white lancet holder cup. (Pushing the lancet into the cup may cock the device; this is okay).

Twist Off the Rounded Top of the Lancet



3. Hold the lancet firmly in place with one hand. Use your other hand to twist off the rounded top of the lancet. Be careful not to touch the exposed needle on the lancet.

Replace the Cap on the Lancing Device



4. Replace the Clear Cap until it snaps or clicks into place. Be careful not to touch the exposed needle on the lancet.

Set the Lancing Level



5. The FreeStyle™ Lancing Device offers five different depth settings. Move the Depth Setting Dial to the desired setting as shown in the Depth Indicator Window. It is recommended that you start at Level 3. Level 1 is the shallowest depth, Level 5 is the deepest.


Cocking the Handle




6. Pull the dark blue Cocking Handle out until it clicks. (You may have already cocked the handle in step 2 above; this is okay).



1. Remove the plastic cover from the top rear of the Visor PDA to expose the slot for the Module.
2. Insert the FreeStyle Tracker Glucose Measurement Module into the slot as shown.
3. “FreeStyle Tracker Diabetes Management System” will appear on the screen, indicating the System is getting ready for a glucose test.

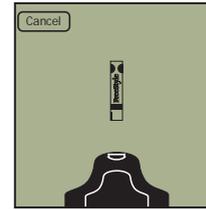


1. This screen will appear after inserting the Module

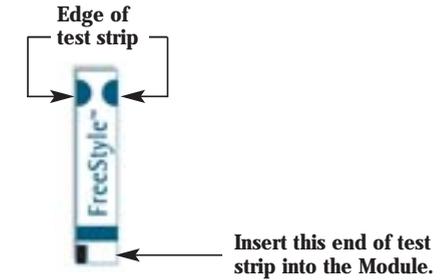

2. Gently insert a new test strip as shown until it cannot go any further, and the Cal Code screen appears.

NOTE: If the Module is already inserted, the glucose measurement function can be started by tapping the  (Measurement) icon in the top, right hand corner on the Logbook or Diary screens (see pages 20-21). Then proceed to the steps above.

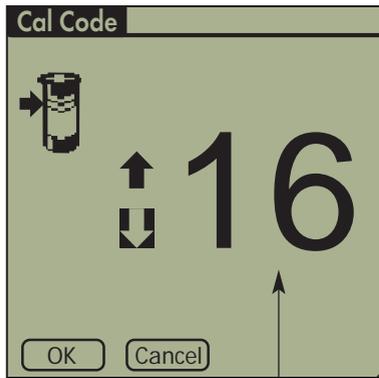

1. The Cal Code screen will appear on the PDA screen after the test strip is inserted.



2. If this Cal Code number does not match the “Code” number in the blue box on the label of the test strip vial, tap the arrows on the screen to adjust to that value.

3. Tap OK when the values match.

NOTE: If you tap Cancel, the logbook screen will appear.

NOTE: Perform a control solution test:

- When you first use the FreeStyle Tracker System, before you take a blood glucose test.
- Each time you open a new vial of test strips.

A control solution test follows the same procedure as a blood glucose test. Substitute a drop of control solution for a drop of blood in performing the test. (For detailed information, go to page 25.)

Important System and Test Strip Information

- The FreeStyle Tracker System has an operating range of 50-95°F (10-35°C).
- Store the test strip package in a cool, dry place below 86°F (30°C). Do not refrigerate or freeze.
- Keep away from direct sunlight or heat.
- Store your test strips in their original vial only; never transfer them to another vial or any other container.
- After removing a test strip from the vial, immediately replace the vial cap and close tightly.
- With clean, dry hands, you may gently touch the test strip anywhere when removing it from the vial or inserting it into the Module.
- Do not use test strips beyond the expiration date printed on the package as this may cause inaccurate results.
- Do not bend, cut, or alter FreeStyle Test Strips in any way.
- Never store individual test strips in your carrying case.

For detailed storage and usage information, refer to the FreeStyle Test Strip package insert.

WARNING

Keep the test strip vial away from children. The cap is a choking hazard. The cap or vial contains drying agents to protect the test strips. Drying agents may be harmful if inhaled or swallowed and may cause skin or eye irritation.

Step 5 – Obtain a Blood Sample

1. Wash your hands and the test site with soap and warm water. Dry both well. Make sure that there is no cream or lotion where you will lance.



2. To bring fresh blood to the surface, rub the area vigorously until you feel it getting warm.



3. Using the FreeStyle Lancing Device (prepared earlier – starting on page 10), press the Clear Cap of the lancing device against the test site.

NOTE: Do not take a sample where veins or moles can be seen. Also, avoid tendons or bony areas.



4. Press the Release Button. **DO NOT LIFT UP.** Continue to hold the lancing device against the test site, gradually increasing the pressure for several seconds.

5. While holding the lancing device on the test site, look through the Clear Cap. The blood sample should be the size of a pin head (example ).
6. Lift the lancing device straight up. Be careful not to touch the blood sample on the test site.

NOTE: It is possible you may experience bruising at the test site. To test at other sites, see the information on page 23.

Step 6 – Apply the Blood Sample

To Apply the Blood Sample to the Test Strip:

CAUTION


- > Use only **ONE** edge of the test strip per test.
- > **DO NOT** “press” the test strip against the test site.
- > **DO NOT** scrape the blood onto the test strip.
- > **DO NOT** apply blood to the flat side of the test strip.
- > **DO NOT** apply blood to test strip when the test strip is out of the Module.
- > **DO NOT** put blood or foreign objects into the strip port.
- > Use the test strip only **ONCE**.



1. Bring the test strip (that has been inserted into the Module) to the blood sample at a slight angle.



2. Gently touch only one edge of the test strip to the blood sample. Do not lift the strip up. Wait until you hear a single beep, or see the Processing screen displayed. This means that enough blood was absorbed.

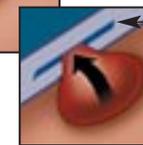

Edge of test strip



Edge of test strip

IMPORTANT

DO NOT apply blood to both edges.

CAUTION

If blood is applied to both edges of the test strip you may get an inaccurate reading. If you have performed a test using both sides of the test strip, redo the blood glucose test with a new test strip and a fresh blood sample.



3. A clock animation on the screen labeled "Processing" indicates that the test is under way.

Inadequate Blood Sample

If the FreeStyle Tracker System has not beeped or the Processing screen does not appear, the sample may be too small. You can continue to fill the test strip for up to sixty seconds. Just be sure you fill the test strip only from the same edge from which you started.

The Result screen will appear when the test has been completed (signaled by 2 beeps).

Result Cal Code = 16 ← Cal Code used for this test

114 ← Test result

2/20/02 7:53 am mg/dL ← Unit of measure

← Tap **New** to start a new test

← Tap **Logbook** to view this result in the Logbook display.

Tap **Details** to add more information to the result

Tap **Control** if this test was done using control solution. A ✓, will be displayed below the Cal Code.

CAUTION

Low or high blood glucose readings (displayed as "Lo" or "Hi") can indicate a potentially serious medical condition. The FreeStyle Tracker System displays results between 20 and 500 mg/dL (1.1 and 27.8 mmol/L).

This concludes the blood glucose test.

- › Discard the used test strip.
- › Test strips may only be used once.
- › To remove the lancet, see page 24.

Viewing Test Results in the Logbook

The Logbook screen contains the following information:

Date → Wed 2/19/02

Time Periods →

	Brkfst	Lunch	Dinner	Bd/Slp
G:	70/170	230*	310/	65/110
I:	4.0	5.0	12.0	16.0*
C:	41	0	80	20

Wed 2/20/02

	Brkfst	Lunch	Dinner	Bd/Slp
G:	60/	210/330	115/	170/130
I:	4.0	7.0	8.0	16.0*
C:	46	100	75	20

G: blood glucose test results. 70/170 in the Breakfast Time Period means 70 was the "Pre Breakfast" result, 170 was the "Post Breakfast" result

An * here indicates this is the earliest of multiple entries - tap it to view the additional entries in the Diary

I: the sum of insulin entries in a Time Period

C: the sum of carbohydrate entries in a Time Period

Tap **New** to make a new entry

Tap the date, ◀▶, or the scroll buttons to change the date displayed

Tap **Today** to bring the lower Logbook chart to today's date

An * here indicates this is the sum of two or more entries - tap it to view the separate entries in the Diary

Tap any entry in the Logbook screen or tap the ☰ icon to go to the Diary screen

Viewing Test Results in the Diary

The Diary screen contains the following information:

Date and time of event → 2/12/02 6:18 a

Event type → <B

Event values → Gluc* 201

Tap an item to highlight it, then tap **Details** to see more information related to the highlighted event


Tap **New** to make a new Diary entry

Tap **Go to** to view events from other dates

Tap these Next / Previous arrows to move forward or back in time in the Diary

This icon indicates there is a comment with this event

Icons and Terms	Explanation
ICONS	
	Go to the Diary screen
	Go to the Logbook screen
	Go to the Statistics screens
	Go to the Line Graph / Pie Chart / Modal Day screens
	Go to the Glucose Measurement screen
TIME PERIODS	
<B	Pre Breakfast
>B	Post Breakfast
<L	Pre Lunch
>L	Post Lunch
<D	Pre Dinner
>D	Post Dinner
Bd	Bedtime
Slp	Sleep
EVENT TYPES	
Gluc	Blood glucose test result (automatic)
Gluc*	Blood glucose test result (manual entry)
Cntl	Control solution test result
Insu	Insulin units taken
Meal	Carbohydrates consumed
Meds	Medication taken
Hlth	State of Health
Exer	Exercise or physical activity performed
Prof	Basal Profile change
Temp	Temporary basal rate change
Note	Additional comments entered
*	In the Logbook, this indicates that there are multiple entries for this item - tap the item to view the other entries

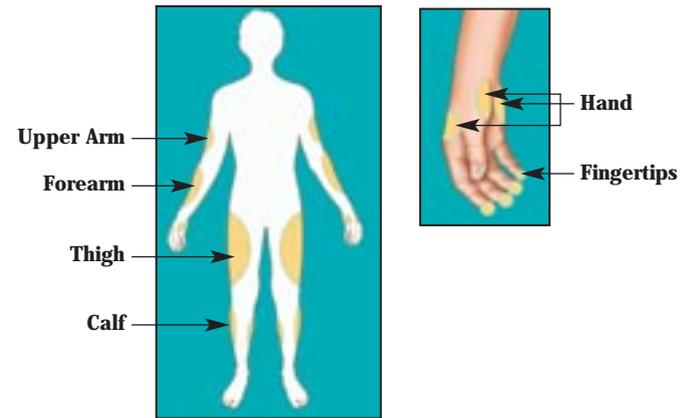

Lancing with the Clear Cap may be completed at the following test sites:

- > Forearm
- > Upper arm
- > Hand
- > Thigh
- > Calf

An example of lancing the forearm is shown on page 16. Testing at other sites is similar, except for the fingers.

IMPORTANT

It is possible you will experience bruising at the listed test sites above. If bruising does occur, you can lance your finger or another site instead. Should you choose to lance your finger, a special dark blue Finger Cap is also provided with your FreeStyle Lancing Device. To test your finger, see the instructions on lancing your finger in the electronic FreeStyle Tracker™ System User's Guide.

WARNING

If you are testing for hypoglycemia (low blood sugar), or if you suffer from hypoglycemia unawareness, we recommend that you test on your fingers.

How to Remove the Lancet

1. When you have finished testing your blood, remove the Clear Cap of the lancing device.
2. Hold the lancet over a sharps container or a puncture-proof container with a lid.



3. Pinch the white clip that holds the lancet until the lancet falls out.

Maintenance

Use isopropyl alcohol or soap and water to wipe the outside of the Lancing Device.

If you wish, remove the Clear Cap. Wash it in warm water and rinse well.

Do not immerse the body of the lancing device in water or any liquid.

When should I perform a Control Solution test?

- When you first get your FreeStyle Tracker™ System, before doing a blood glucose test.
- When you first begin using a new vial of test strips.
- When you suspect that your FreeStyle Tracker System or Test Strips are not working properly.
- When you think your test results are not accurate, or if your test results are not consistent with how you feel.
- When your FreeStyle Tracker System is dropped, damaged, or exposed to liquids.
- When your healthcare professional advises you to do so.

Precautions when using Control Solution

- Use only FreeStyle brand control solution with your FreeStyle Tracker System. Use of other brands of test strips or control solution with the FreeStyle Tracker System may produce inaccurate results.
- Replace the cap on the control solution bottle immediately after using.
- Results from FreeStyle Control Solution tests do not reflect your blood glucose levels.
- FreeStyle Control Solution is good for 3 months after opening the bottle, or until the expiration date printed on the label, whichever occurs first.
- Do not use control solution past the expiration date or you may get inaccurate results.
- Do not add water or any liquid to control solution.
- Control solution tests are specified to be accurate only when used between 59-104°F (15-40°C).

How do I know if my Control Solution has expired?

On the date you open a new control solution vial, count forward 3 months. This is your discard date. Write this date on the side of the control solution vial.

Example: You open the control solution on January 15th. Count forward 3 months to April 15th. This is the discard date that you record on the vial.

IMPORTANT

The control solution range printed on the test strip vial label is a target range for control solution only. It is not a target range for your blood glucose level.

If the control solution results are out of this range, REPEAT THE TEST.

Out of range test results may be caused by:

- > Expired control solution
- > Expired or bad test strip
- > Error in performing the test
- > Diluted or watered-down control solution
- > Code on the test strip vial does not match the Cal Code set on the FreeStyle Tracker PDA screen
- > System malfunction
- > The control solution test being done outside the temperature range of 59-104°F (15-40°C)

If your control solution test results continue to fall outside the range printed on the test strip vial, the FreeStyle Tracker System may not be working properly.


- > IF SO, DO NOT USE THE FREESTYLE TRACKER SYSTEM TO TEST YOUR BLOOD.
- > CALL THERASENSE CUSTOMER SERVICE AT 1-888-522-5226.

System Does Not Work During a Glucose Test	
PROBABLE CAUSE	WHAT TO DO
Batteries are dead.	Replace the batteries.
Batteries are installed incorrectly, or there are no batteries in the system.	Check that both batteries are installed correctly. (Refer to PDA instructions.)
Test strip is inserted upside down, wrong end in or incompletely inserted into the FreeStyle Tracker Module.	Insert the test strip correctly with the name “FreeStyle” facing up, and the edge of the test strip clearly visible. (See page 13.)
Defective Module or test strips.	Call TheraSense Customer Service.
Blood or foreign objects put into test strip port.	Call TheraSense Customer Service.

Test Does Not Start After Applying The Blood Sample	
PROBABLE CAUSE	WHAT TO DO
Blood sample too small.	Repeat the test using a new test strip and a larger blood sample.
Defective test strip.	Repeat the test with a new test strip.
Sample applied after automatic shutoff (2 minutes).	Repeat the test using a new test strip. Wait until you see the “Blood Drop and Test Strip” animation on the PDA screen before you apply blood.
Defective system or test strips.	Call TheraSense Customer Service.

The following Error Messages may be displayed if there is a problem with the System or a measurement result. Possible causes for the error are listed, followed by actions you can take to resolve the problem.

Error 0

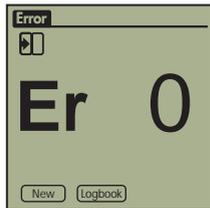


Possible Cause

- > Problem with the Module or
- > Problem with the PDA or
- > Module not properly inserted into the PDA

What You Should Do:

1. Unplug the Module and plug it back into the PDA. Then retest using control solution and a new test strip. If the control solution test works, retest using blood and a new test strip.
2. If the retest using control solution does not work, or if the error persists, call TheraSense Customer Service.

Error 1

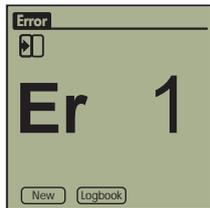


Possible Cause

- > Sample too small or
- > Problem with the test strip or
- > Problem with the system or
- > VERY low blood glucose (less than 20 mg/dL).

What You Should Do:

1. If you have symptoms such as weakness, sweating, nervousness, headache, or confusion, follow your doctor's recommendation for treating hypoglycemia (low blood glucose).
2. Retest using control solution and a new test strip. If the control solution test works, retest using blood and a new test strip.
3. If the retest using control solution does not work, or the error persists, call TheraSense Customer Service.

Error 2

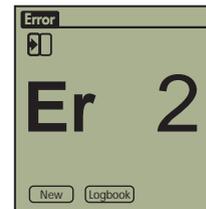


Possible Cause

- > Problem with the test strip or
- > Problem with the system or
- > VERY high glucose (above 500 mg/dL).
- > HI control solution applied when temperature is too cold (condensation). (This applies only to control solution that is labeled as HI).

What You Should Do:

1. If you have symptoms such as thirst, fatigue, excess urination, or blurry vision, follow your doctor's recommendation for treating hyperglycemia (high blood glucose).
2. Retest using control solution and a new test strip. If the control solution test works, retest using blood and a new test strip.
3. If the retest using control solution does not work, or the error persists, call TheraSense Customer Service.

Error 3

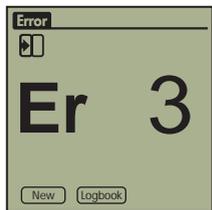


Possible Cause

- > Incorrect test procedure, e.g. putting blood on the test strip before inserting the test strip in the Module, blood applied before the “Blood Drop and Test Strip” animation or
- > Problem with the test strip or
- > Problem with the System.

What You Should Do:

1. Be sure you see the “Blood Drop and Test Strip” animation on the PDA screen before you apply blood or control solution.
2. Retest using control solution and a new test strip. If the control solution test works, retest using blood and a new test strip.
3. If the retest using control solution does not work, or if the error persists, call TheraSense Customer Service.

Error 4

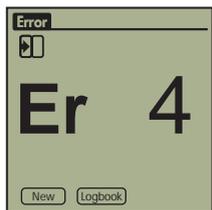

Possible Cause

- > Problem with the test strip or
- > System malfunction.

What You Should Do:


1. Retest using control solution and a new test strip. If the control solution test works, retest using blood and a new test strip.

2. If the retest using control solution does not work, or if the error persists, call TheraSense Customer Service.

Error 5

Error 5 is not used with the FreeStyle Tracker System.

Error 6

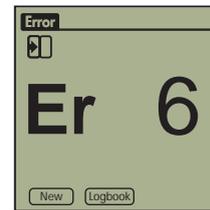

Possible Cause

- > A button on the system was pushed when applying blood sample or control solution

What You Should Do

1. Retest using a new test strip. Do not touch any buttons when applying blood or control solution to the test strip.
2. If the error persists, call TheraSense Customer Service.

Avoid getting dirt, dust, blood, control solution, water, or any liquid inside the FreeStyle Tracker Module through the test strip port and the Module connector port.

Store the System in the carrying case after each use.

Specifications	
FREESTYLE TRACKER SYSTEM TESTING SPECIFICATIONS	
Result Range	20 to 500 mg/dL (1.1 - 27.8 mmol/L)
Calibration	Plasma equivalent
Sample	Whole blood, capillary
Sample Size	300 nanoliters (0.3 microliters)
Test Time	Average of 15 seconds
Assay Method	Coulometric Electrochemical Sensor
Measurement Units	Either mg/dL or mmol/L
Module Memory Backup	250 most recent blood glucose and control solution tests
PDA Logbook Memory	2,500 glucose tests and other entered events
Automatic Shut-off	2 minutes and user adjustable
PDA Power Source	2 replaceable AAA batteries, or rechargeable PDA battery
AAA Battery Life	2 - 4 weeks
Operating Temperature	50 - 95 °F (10 - 35 °C)
Operating Relative Humidity	5-90%
Hematocrit	0-60%



1360 South Loop Road
Alameda, CA 94502

www.freestyletracker.com
Toll Free 1-888-522-5226